
# Invia® Wound Therapy



## WOUND DRESSING SETS

EN Instructions for Use – Internal Drain/External Interface Application Instructions


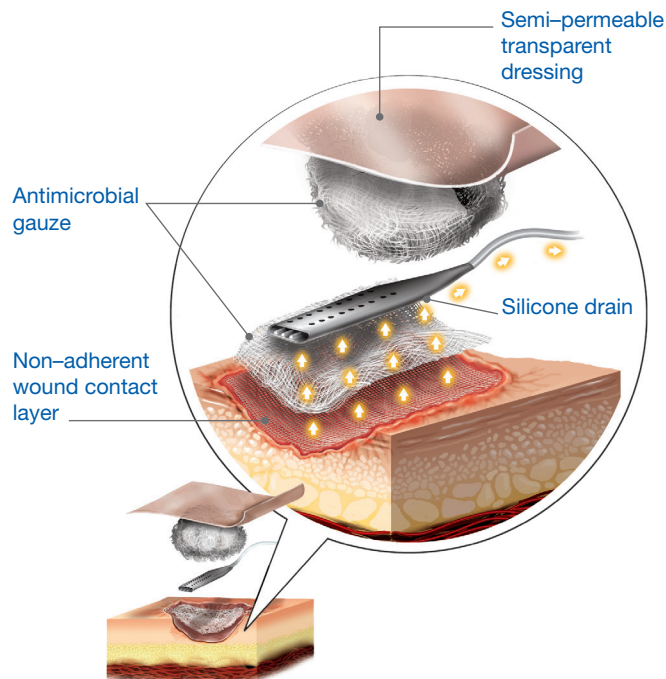
Please keep in mind that each wound is unique and must be assessed by a qualified medical professional who has read and adheres to the Instructions for Use.

1. Inspect the dressing set for completeness and ensure there is no damage to packaging or individual items within the kit.
2. Cleanse wound and periwound area according to institution guidelines or physicians orders. Irrigate the wound thoroughly with saline and gently wipe away any excess fluid with sterile antimicrobial gauze.
3. Prep the skin up to 2 inches around the wound bed using the Cavilon™ skin sealant. Allow to dry.
4. Cut non-adherent wound contact layer to mimic the size and shape of the wound bed. Cover as much of the wound as possible with the non-adherent dressing. (See Illustration A)

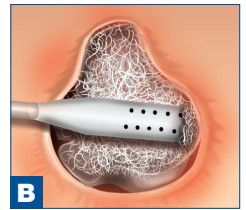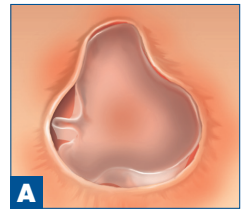Flat drain pictured, other styles available.

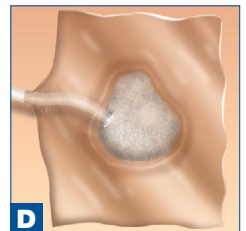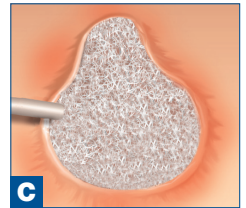Note: If you are treating two wounds simultaneously, a Y-connector is available from Medela.

### Internal Drain Instructions:

5. Trim the drain (See Illustration B) so that it fits within the wound bed, then either:
  - a. Wrap the drain in antimicrobial gauze and place it in the wound on top of the non-adherent contact layer; or
  - b. Saturate antimicrobial gauze with saline place it in the wound bed top of the non-adherent contact layer and then place the drain directly on top of the gauze.
6. Saturate additional antimicrobial gauze with saline, open and fluff into the wound to completely cover the drain and fill the wound to skin level. (See Illustration C)
7. Cut transparent dressing so that it will extend at least 2 inches beyond the wound edge on all sides.
8. Create an airtight seal by applying the transparent dressing starting with the side opposite the tubing. Adhere the transparent dressing firmly to the entire wound site. When you reach the tubing, you may lift it up and squeeze the dressing around and under the tube to create an airtight seal.
9. Press the transparent dressing down firmly around the entire wound to ensure the wound environment is sealed. (See Illustration D)
10. Strip paste may be used to comfortably affix tubing and to promote a good seal in difficult to dress areas, while protecting the integrity of the skin. (See Illustration E)



11. Connect the end of the selected drain to the double lumen tubing conical connector assuring that the all clamps are unclamped and turn on the pump with desired physician ordered settings.

12. The dressing will compress when pump starts and dressing will contract noticeably and feel firm to touch. If the dressing does not contract, the system may not be sealed. Inspect the dressing and reseal if necessary.

13. Latex-free, hypoallergenic surgical tape may be used to window pane transparent dressing and affix drainage tubing as necessary.

14. Dressing changes are recommended 2 to 3 times weekly or as directed by a physician.

15. Routinely check that the vacuum level is set as prescribed and that the dressing is properly sealed.
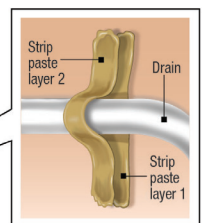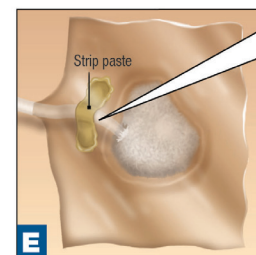
## Instructions for Use — External Interface Application Instructions

1. Inspect the dressing set for completeness and ensure there is no damage to packaging or individual items within the kit.
2. Cleanse wound and periwound area according to institution guidelines or physicians orders. Irrigate the wound thoroughly with saline and gently wipe away any excess fluid with sterile antimicrobial gauze.
3. Prep the skin up to 2 inches around the wound bed using the Cavilon™ skin sealant. Allow to dry.
4. Cut non-adherent wound contact layer to mimic the size and shape of the wound bed. Cover as much of the wound as possible with the non-adherent dressing. (See Illustration F)
5. Saturate antimicrobial gauze with saline, open and fluff into the wound to completely fill the wound to skin level.
6. Cut transparent dressing so that it will extend at least 2 inches beyond the wound on all sides, and grasp transparent film by either tabbed end and remove backing.
7. Holding a tabbed end apply the transparent dressing starting on one side of the wound, apply to the skin smoothing and sealing as it is applied, remove remaining tabs. Place gloved hands over the transparent film for a few seconds, press down firmly to assist adherence to the entire wound site and ensure a good seal.
8. Select an appropriate location over the dressed wound bed for the ESI to be applied and pinch the transparent film over the area lifting slightly, and cut a hole in the pinched area approximately 1 cm in width. The hole should not exceed the size of the ESI to prevent a leak in the system. (Illustration F)
9. Peel the backing off of the ESI and center the ESI over the previously cut hole in the dressed wound and apply, pressing firmly for adherence. Should additional tape be desired it may be applied in a window pane fashion over the perimeter of the transparent film and the perimeter of the ESI. (Illustrations G & H)
10. Connect the end of the selected drain to the double lumen tubing conical connector assuring that the all clamps are unclamped and turn on the pump with desired physician ordered settings.
11. The dressing will compress when pump starts and dressing will contract noticeably and feel firm to touch. If the dressing does not contract, the system may not be sealed. Inspect the dressing and reseal if necessary.
12. Latex-free, hypoallergenic surgical tape may be used to window pane transparent dressing and affix drain or drainage tubing as necessary.
13. Dressing changes are recommended 2 to 3 times weekly or as directed by a physician.
14. Routinely check that the vacuum level is set as prescribed and that the dressing is properly sealed.

**Invia® Wound Therapy Dressing Sets can be used with Invia Motion™ and Invia Liberty™ Negative Pressure Wound Therapy Systems. Please consult the Instructions for Use for the specific system that is being used with the dressing set to understand the applicable indications for use, contraindications, warnings, and cautions.**

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